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

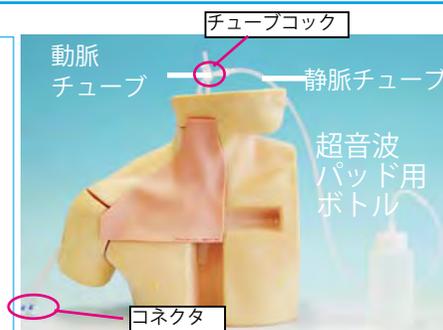
Instruction Manual
 Ultrasound Puncture Pad for M93U CVC Simulator

SET UP

- 1) Place the pad into the torso.
- 2) Put the torso upright.
- 3) Fill the artery tube
- 4) Fill the vein tube
- 5) Lay the simulator and start training.

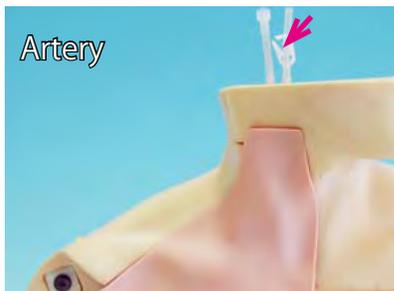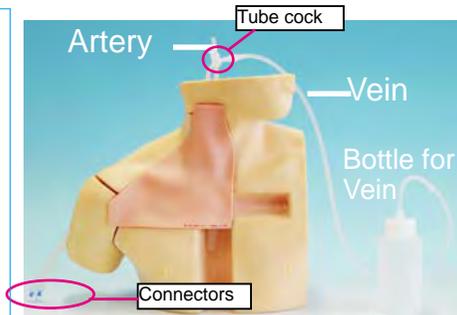
First, fill the artery tube.

Put the torso upright.

This is to avoid bubbles being caught in the tube, as well as to prevent the water flows out from its open end.

Open the tube cock at the head end of the artery tube

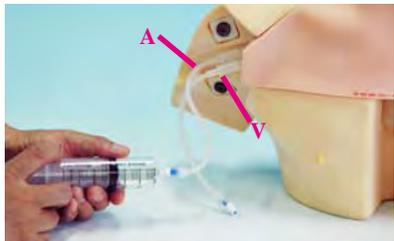Connect the water filled syringe to the connector at the lower end of the artery tube. (Screw in the tip of the syringe to the connector)

Advance the piston of the syringe slowly until the water surface reaches at the height of the tube cock.

When the tube has been filled by water, close the tube cock.

Take off the syringe from the tube.

Then, Fill the vein

Connect the tip of the tube from the bottle to the head end of the artery tube.

Connect the water filled syringe to the connector at the lower end of the vein tube.

Advance the piston of the syringe slowly until the water flows out into the bottle.

Keep the bottle connected while the simulator is in use.

Lay the simulator and start your training session.

■ Main Office and Factories (World Wide)

KYOTO KAGAKU co.,LTD

WEB • www.kyotogagaku.com
 E-MAIL • rw-kyoto@kyotokagaku.co.jp

15 Kitanechoya-cho Fushimi-ku Kyoto 612-8388, Japan
 TEL: +81-75-605-2510
 FAX: +81-75-605-2519

■ Americas

KKAmerica
 Kyoto Kagaku America Inc.

WEB • www.kkameric-inc.com
 E-MAIL • info@kkamrica-inc.com

3109 Lomita Boulevard, Torrance, CA 90505-5108, USA
 TEL: +1-310-325-8860
 (Toll-free in North America: 877-648-8195)
 FAX: +1-310-325-8867