



Great addition as a demonstration tool for airway management

Compact table-top design, anatomical accuracy and a variety of DAM settings. Ideal for demonstration of skills and devices.

Features

1. Anatomically correct airway
2. The incisors are removable when excessive force is applied.
3. Successful tube tip placement can be confirmed by indicators.

Training Skills

- Airway opening techniques (head tilt, jaw thrust)
- Bag-Valve-Mask ventilation
- Pre-intubation airway assessment
- "Sniffing position"
- Pressurization of external larynx to improve the laryngeal view
- Intraoral/Intranasal Intubation
- Use of oropharyngeal airway (OPA)
- Use of nasopharyngeal airway (NPA)
- Use of laryngeal mask airway
- Use of video laryngeal scope
- Confirmation of successful ventilation by indicators.
- Feedback of incorrect procedures including esophagus intubation and unilateral intubation
- Securing the tube in place with tape or Thomas™ endotracheal tube holder

Recommended Devices:

- (Macintosh laryngoscope) blade size 3,4
- (Endotracheal tube) internal diameter 7.0mm, 7.5mm
- (Laryngeal mask) Air-Q: size 3.5
LMA Supreme: size 4
I-gel: size 4

Specifications :

Set Includes:	Manikin Size:
1 manikin	W45 x D20 x H25 cm
3 upper incisors	
1 lubricant	Replacement Parts:
1 syringe	(11393-020) 1 face mask
1 carrying bag	(11392-080) 2 replacement tongues II
1 instruction manual	(11392-010) 10 upper incisors



Realistic trainings in insertion and management of a fiberscope


Anatomically correct trachea and bronchi provides realistic view through a scope, allowing recognition of bifurcations.

Features

1. Anatomically correct trachea and bronchi that include trachea bifurcation for segmental bronchi.

Training Skills

- Insertion and management of a fiberscope
- Recognition of bifurcations by a fiberscope
- Airway opening techniques (head tilt, jaw thrust)
- "Sniffing position"

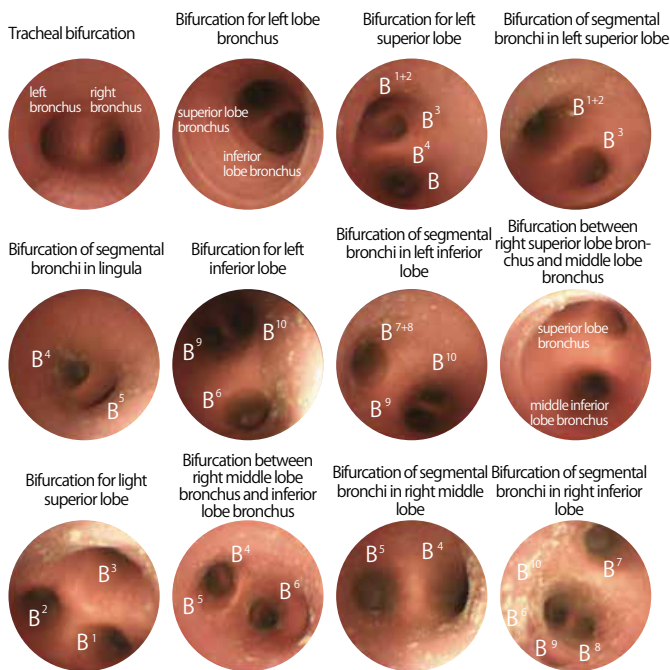

Specifications :

Set Includes:	Manikin Size:
1 manikin	W45 x D20 x H25 cm
3 upper incisors	
1 lubricant	Replacement Parts:
1 syringe	(11393-020) 1 face mask
1 carrying bag	(11392-080) 2 replacement tongues II
1 instruction manual	(11392-010) 10 upper incisors