



**External genitalia made of the new material and anatomically correct vaginal canal provides true-to-life hands-on training in cervical examination to determine Bishop score from the final stage of pregnancy to the first stage of the labor.**

### Features

1. Genital area made of innovative material that excels in elasticity, durability, and resilience, offers repetitive training.
2. Five Interchangeable cervical dilation inserts can be set by a single action.
3. Anatomically correct landmarks including the ischial spines and the pubic symphysis
4. Anterior/posterior fontanelles to determine the fetal position

### Training Skills


- Cervical examination
  - ▶ Determination of Bishop score (cervix score) using landmarks of the ischial spines and the pubic symphysis
  - ▶ Assessment of cervical dilation
  - ▶ Procedures of urinary catheterization
  - ▶ Procedures of cervical cytology
- Speculum examination

### Key Features

**Eternal genitalia that excels in elasticity, durability and resilience**



**Ischial spines and pubic symphysis to determine fetal station**

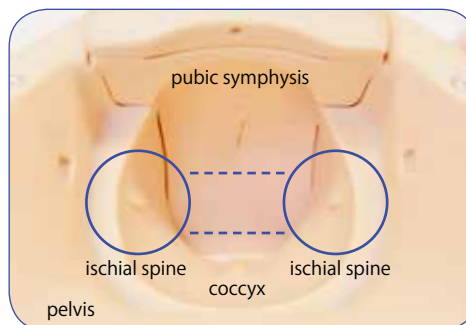

Station zero, the line that connect ischial spines can be recognized.



Cervical dilation insert can be attached by a single action.



Descent of the fetus can be set by the lever on the holder. (from -3 to +3)

### Specifications :

Set Includes:		Manikin Size:	Replacement Parts:
1 maternal body	1 holder for the cervical dilation inserts	W74 x D44 x H33.5 cm 10.3kg	(11415-010) 1 set of cervical dilation inserts (5 variations)
1 genitalia unit	3 lubricant (150 ml)		(11415-020) 3 lubricant (150 ml)
1 vagina unit	1 talcum powder	Materials	(11415-030)
1 set of cervical dilation inserts (5 variations)	1 instruction manual	Soft resin, hard plastic	1 genitalia unit (for pelvic examination and delivery assistance)