


# Prostate & Rectal Simulator

MW20 11399-000

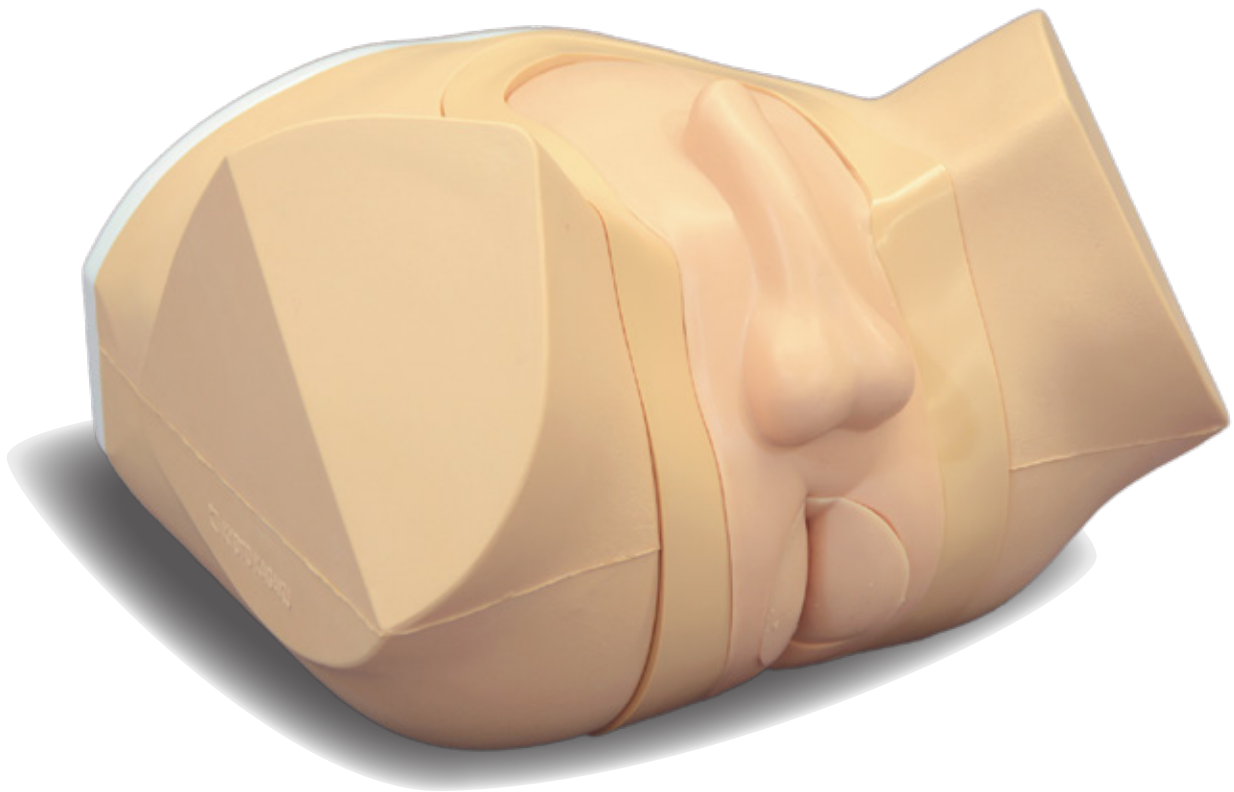

SET INCLUDES		SIZE
3 prostate units	1 lubricant gel	W34 x D33 x H24 cm
: 9 prostates	1 pelvic body	13.4W x 13D x 9.45H inch
4 rectum units	1 storage case	

Specifications are subject to change.

 **KYOTO KAGAKU co.,LTD**

Worldwide Inquiries & Ordering

<http://www.kyotokagaku.com>

E-mail: [rw-kyoto@kyotokagaku.co.jp](mailto:rw-kyoto@kyotokagaku.co.jp)

**Kyotokagaku Head Office and Factories**

15 Kitane-koya-cho, Fushimi-ku, Kyoto, 612-8388, JAPAN

Tel: +81-75-605-2510 Fax: +81-75-605-2519

# Prostate & Rectal Simulator


The Prostate & Rectal Simulator include 9 types of prostates and 4 types of rectums that are easily changeable and can be positioned in 3 different ways.

## FEATURES

- Nine types of prostate simulate different cases
- Four rectum units
- Training with multiple prostate and rectal cases is concurrent and time efficient
- Lateral, supine and prone positioning
- 3 revolving units with 9 prostates can easily be changed.

## TRAINING SKILLS

- Digital-examination of the prostate and rectum
- Insertion and the use of the anal-scope and the proctoscope

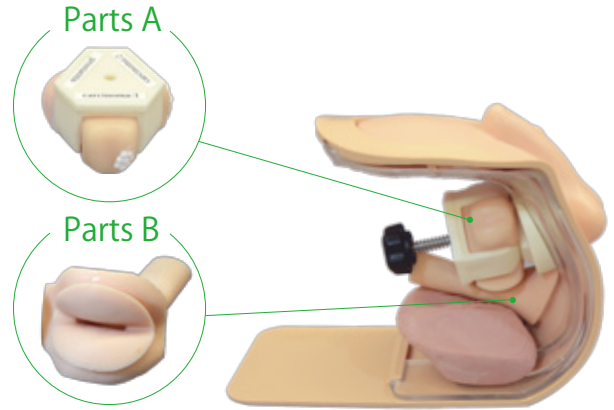


## Parts A Nine Types of Prostates



- |                  |                |                |
|------------------|----------------|----------------|
| #1 normal        | #4 prostatitis | #7 carcinoma 3 |
| #2 enlargement 1 | #5 carcinoma 1 | #8 carcinoma 4 |
| #3 enlargement 2 | #6 carcinoma 2 | #9 carcinoma 5 |

The hardness of each prostate is different. When in use the prostate is hidden and cannot be seen by the trainee. 3 revolving units and 9 prostates can easily be changed.

## Parts B Four Rectum Units

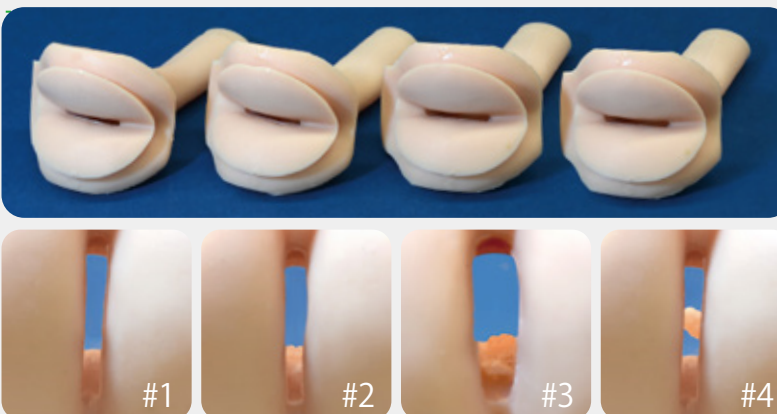

- #1 normal
- #2 small carcinoma
- #3 large carcinoma
- #4 polyp combined with small carcinoma